

# Product datasheet

## Quadroflex Luxe

*A 4-section table; as versatile as you. With hinges beneath your patient's knee, hip and shoulder joints, the Quadroflex provides ideal support. This anatomically based split makes almost every position possible. Ideal in case of moving your patient from laying to sitting position. The patient does not have to rotate his/her feet on the head section to take a sitting position which raises comfort and hygiene.*

*Positioning with the Quadroflex is easy thanks to the 3 powerful motors.*

### Standard features

- Four-section table with adjustable head-, back-, sitting-, and leg section.
- A total of 3 motors control height, back rest and sitting position.
- Gymna cushion heating with 2 temperature settings and hand control.
- Ergomax Face Support for optimum patient comfort while lying in prone position (one supplied as standard).
- Arm supports in the head section.
- Optional fold-down supports in the back section.
- Electrical height adjustment, operated by peripheral foot switch from all 4 sides of the table.
- Gymna Stability Profile (GSP).
- Cushions with duo-padding, with a high-density lower layer.
- The Luxe cushions are finished with tough saddle stitching.
- Central wheel-lifting mechanism with 4 casters.
- Non slip feet.
- Cushions in a choice of 11 colours.
- Prepared for shoulder traction supports.

### Ordering data

125 094 Quadroflex Luxe incl. peripheral switch and cushion heating

### Specifications

Length x width:	± 202 x 67cm
Height adjustment:	± 45 – 95 cm
Weight:	± 112 kg
Max. lifting capacity:	200 kg
Electric motor:	6.000 N
Cushion temperature:	35 or 38 °C
Power supply:	100 – 240 VAC, 50-60 Hz
Current consumption:	3,0 A
Insulation class:	II
Moister resistance:	IP 21
Classification:	I
Conformity:	IEC/EN 60601-1 IEC/EN 60601-2 CAN/CSA C22.2 No. 601.1 MDD 93/42/EEC
Guarantee:	3 years (upholstery 1 year)



### Optional

125 622 Fold-down sections

### Accessories

150 414	Ergomax Face Support
122 405	Paper roll holder, for undercarriage
337 018	Paper rolls (per 12 pieces)
107 819	Traction support for shoulders
121 401	Gynaecological leg support
333 245	Traction device AutoTrac 460
121 434	Attachment frame for AutoTrac 450/460
332 915	Foot end protector for cushion
xxx xxx	Terry-cloth cover, stretch material (available in 9 colours 'Style')



Gymna Stability Profile



Gymna cushion heating inside


## Quadroflex Luxe



**NEW:** Ergomax Face Support

*Quadroflex with optional fold down section*

### Functionality/dimensions

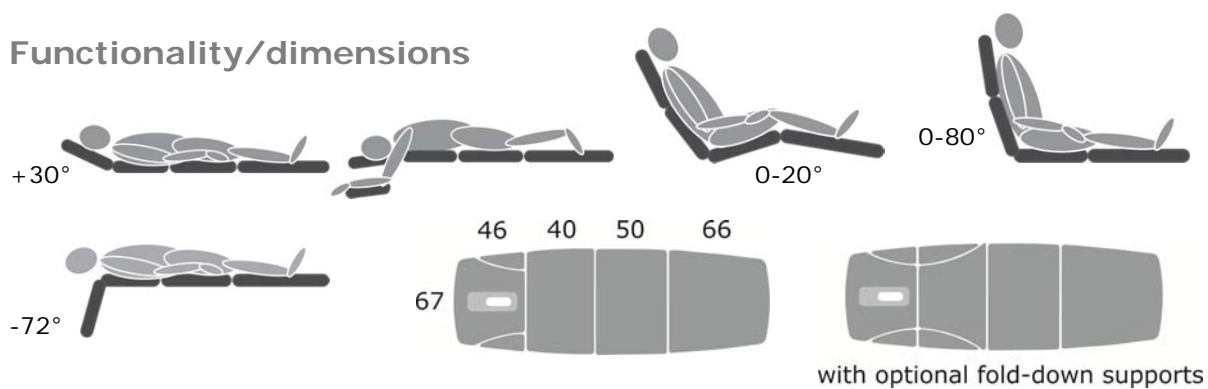

### Cushion colours

#### Style



#### Fashion

